

Princess 360 Water Birth Pool

With a footprint of just 1.76m² the freestanding Princess 360 is ideal for rooms where space is at a premium.


Mothers and midwives benefit from complete freedom of movement both in and around the pool and will appreciate the comfort and support they receive.

The Princess 360 is perfectly formed, surprisingly spacious and delivers the same performance that our larger pools do.



Specifications

Length:	1600mm
Width:	1100mm
Height:	770mm
Depth:	640mm
Weight (empty):	70kg
Capacity:	450 litres
Material:	Ficore® composite
Guarantee:	Lifetime

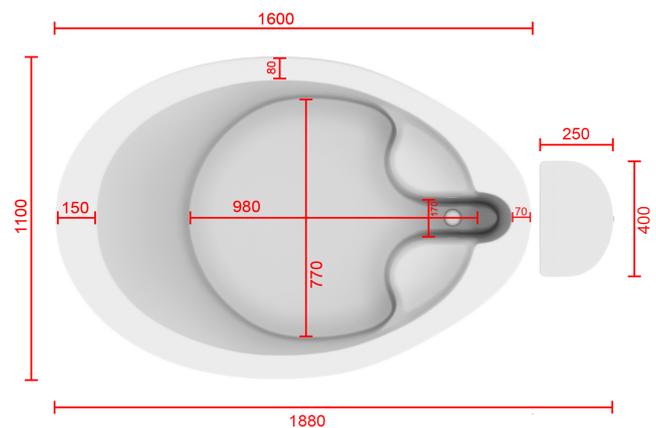

Princess 360 Water Birth Pool

ACTIVE
BIRTHPOOLS

The Labour Support Seat encourages the pelvis to open and provides support in upright postures

The Water Column separates the plumbing from the pool to guarantee highest level of hygiene and water safety

The keyhole shaped channel in the rim protects mothers from hitting their head on the spout

Integral LED Lighting and Bluetooth Sound allow mothers to personalise the ambiance of the room

The step unit and grab rail help make it easy for mothers to safely get in and out of the pool

The wrap-around rim provides support for mothers in the full range of positions natural to labour and birth

The highly polished finish is less slippery, super smooth, tactile and warm to the touch

360 degrees of movement around and in the pool for mothers and midwives

The Lumbar Support provides comfort and support for mothers in upright postures like squatting

High density and superior strength minimise heat loss and make our pools more durable and economical

Though compact, the Princess Pool is surprisingly spacious for mothers to move freely in comfort

Designed to safely facilitate emergency evacuations with or without a portable hoist

Easy to clean, resistant to bacteria and repairable on site, like new if damaged

There's plenty of space under the rim for midwives to get close to the mother in comfort

