


TECHNICAL SPECIFICATIONS

ACTIVE
BIRTHPOOLS
www.activebirthpools.com

Active II/360 Birth Pool

Material	Ficore composite
Capacity	560 liter / 145 gal (4" below rim)
Height	770mm / 30.25"
Width	1650mm / 65"
Length	2260mm / 89"
Water Depth	640mm / 25" (4" below rim)
Weight (empty)	100kg / 220 lbs
Shipping Crate	205 x 181 x 97cm / 81" x 71" x 38"
Shipping Weight	220kg / 485 lbs
Birth Stool	2 - 1 @ water level, 1 submerged
Handles	2 - Chrome, permanent (Electrical grounding available upon request)
Drain Mechanism	40mm / 1 1/2" Mechanical pop-up—brass with cable mechanism, rim mounted control, removable and height adjustable chrome stopper
Waste Plumbing	40mm / 1 1/2" PVC - Shallow P-trap included - thru-wall (or thru-floor connection - omit supplied p-trap)
Access Panel	2 - right and left rear of pool - 10mm hex key lock
Chromotherapy Light (optional)	1 - 9 light, 11 color LED: White, light blue, dark blue, purple, magenta, red, pink, orange, yellow, light green, dark green
Transformer (optional)	120VAC to 12VDC 0.5A - CSA and UL Certified - NEMA 5-15 plug
Light Controls (optional)	Rim mounted - chrome, sealed, four button - Mounted on right side unless specified
Step Stool (optional)	18" x 12" x 6", 8.5 lbs, Ficore composite

www.activebirthpools.com

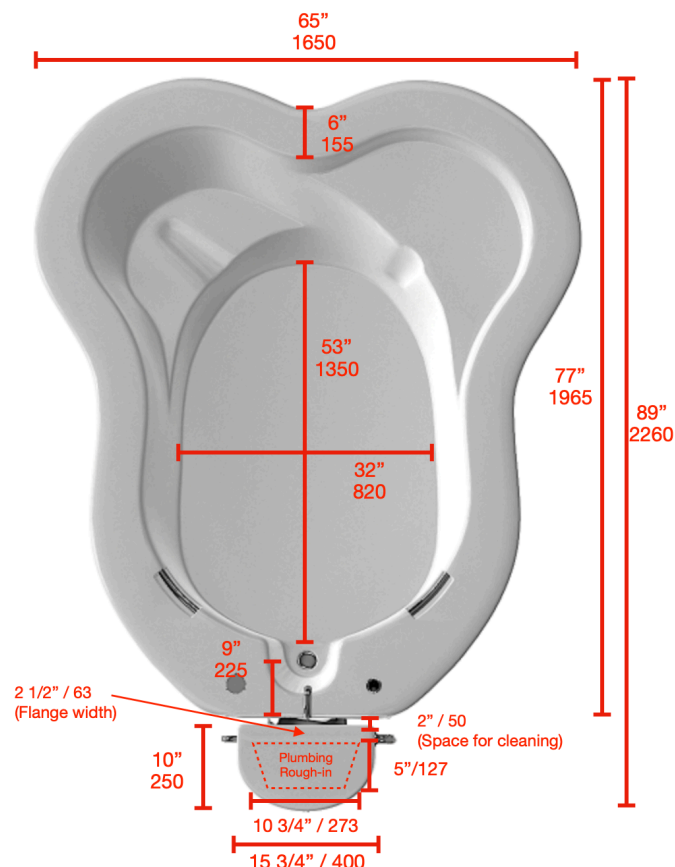

NOTE: Dimensions are in millimeters if not indicated otherwise.

Active II/360 Water Birth Pool

ACTIVE
BIRTHPOOLS

The bull-nose shaped rim provides firm hand-hold for support in upright positions

Multi-colour LED lighting and Bluetooth sound system allow mothers to personalise the ambiance of the room

The sturdy column provides support for upright birth positions and a plinth for midwives to lean on

The keyhole shape channel in the back rim protects mothers from hitting their head on the spout

The recessed hand grips provide the ideal support for upright positions such as squatting

The wide rim, grab rail and step unit make it simple and safe for mothers to get in and out of the pool

Designed to safeguard mothers and midwives from risk of injury and cross infection

The rounded, bull nose-shaped rim provides dynamic support for upright positions

Specialised material makes our pools less slippery, tactile and warm to the touch

This large, luxurious pool gives mothers and midwives complete freedom of movement both in and around the pool

The extra hard, dense material is highly resistant to damage, bacteria and disinfectants

Midwives sit naturally, in comfort with their legs under the pool and arms resting on the extra-wide rim

The concave, cutaway surround enables portable hoists to be utilised for emergency evacuations

The innovative Safety Seat provides support for mothers and gives midwives safe, practical options for monitoring, and emergency evacuations

Mothers have optimum space, depth and support to gain maximum benefit from water immersion

The unique Labour Support Seat is tilted forward to encourage the pelvis to spread and the birth canal to open.

