


TECHNICAL SPECIFICATIONS

ACTIVE
BIRTHPOOLS
www.activebirthpools.com

Venus 360 Birth Pool

Material	Ficore composite
Capacity	510 Liter / 135 Gallons (4" below rim)
Height	770 mm / 30"
Width	1350 mm / 53"
Length	2345 mm / 92"
Water Depth	650 mm / 25" (4" below rim)
Weight (empty)	90kg / 198 lbs
Shipping Crate	77" x 56" x 35"
Shipping Weight	407 lbs
Birth Stool	1 - submerged
Handles	2 - Chrome, permanent (Electrical grounding available upon request)
Drain Mechanism	40mm / 1 1/2" Mechanical Pop-up—brass with cable mechanism, rim mounted control, removable and height adjustable chrome stopper
Waste Plumbing	40mm / 1 1/2" PVC - thru-floor connection only
Access Panel	1 - rear of pool - 10mm hex key lock
Chromotherapy Light	1 - 9 light, 11 color LED: White, light blue, dark blue, purple, magenta, red, pink, orange, yellow, light green, dark green
Transformer	120VAC to 12VDC 0.5A - CSA and UL Certified - NEMA 5-15 plug
Light Controls	Rim mounted - chrome, sealed, four button - Mounted on right side unless specified
Step Stool	18" x 12" x 6", 8.5 lbs, Ficore composite

www.activebirthpools.com

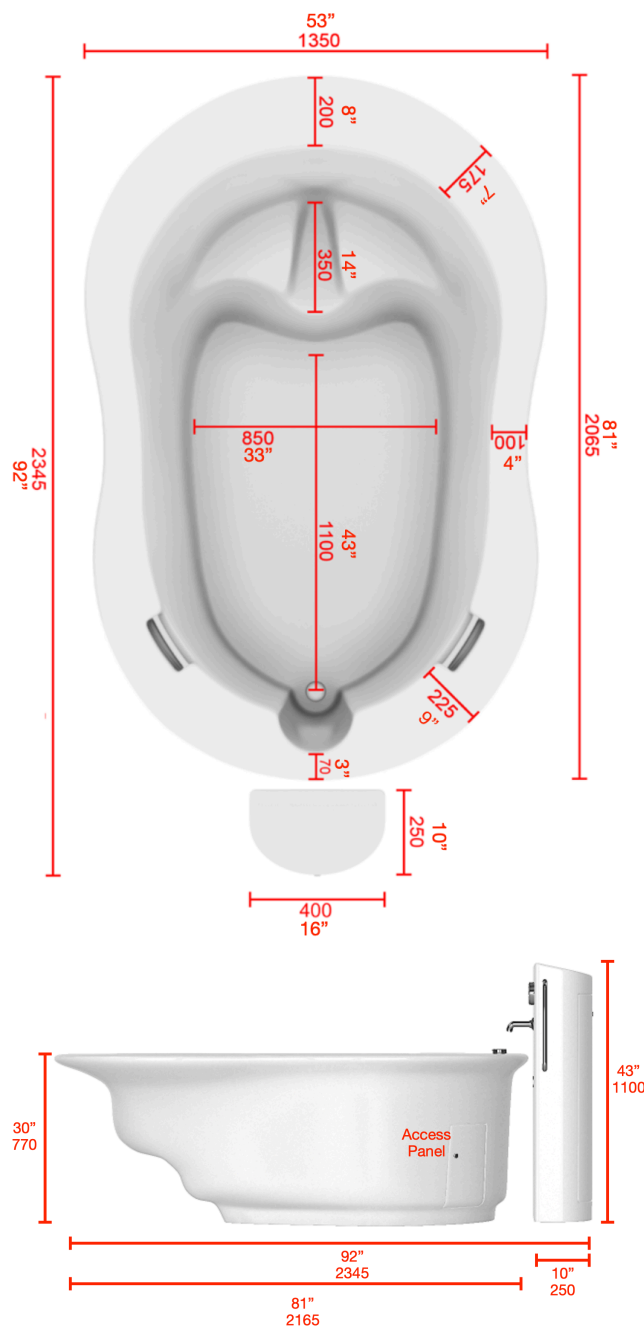

- Water Column Sold Separately -

NOTE: Dimensions are in millimeters if not indicated otherwise.

Venus 360 Water Birth Pool

ACTIVE
BIRTHPOOLS

Multi-colour lighting and Blue-tooth sound enable mothers to personalise the ambience of the room

The Water Column separates the plumbing from the pool to maximise hygiene and water safety

The 'keyhole' shaped channel protects the mother from hitting her head on the spout

The recessed hand grips are bonded into fabric of the pool to deny bacteria to a place to grow

The extra-wide rim, grab rail and step unit make it easy for mothers get in and out of the pool

The combination of double wall construction and ultra-dense material reduces heat loss to .7 degrees per hour

Provides freedom of movement for mothers and midwives to move naturally both in and around the pool

Seamless one-piece construction in Ficore composite makes the pool more durable and easier to clean

The Labour support Seat encourages optimal positioning and gives mothers a place to rest and welcome their baby

Mothers use the extra-wide rim to rest on and to safely enter and leave the pool

Easier to clean and maintain and can be easily repaired on site like new if damaged

The concave, cut-away surround makes it easier and more comfortable for midwives to be close to the mother

Optimal space and depth for mothers to move and find most beneficial positions

The bull-nose shaped rim provides mothers with support in the wide range of postures natural to labour and birth

